

re|Fill

face & skin lifting

feel the difference

re|Fill




Face Lifting – Experience your new identity

Comfortable, Effective, Automatic

The new generation of lifting treatments. reFill ensures a homogen result by automatic injection with a 9-Pin needle adapter.

From small wrinkles to bigger wrinkles reFill will support the doctor's needs.

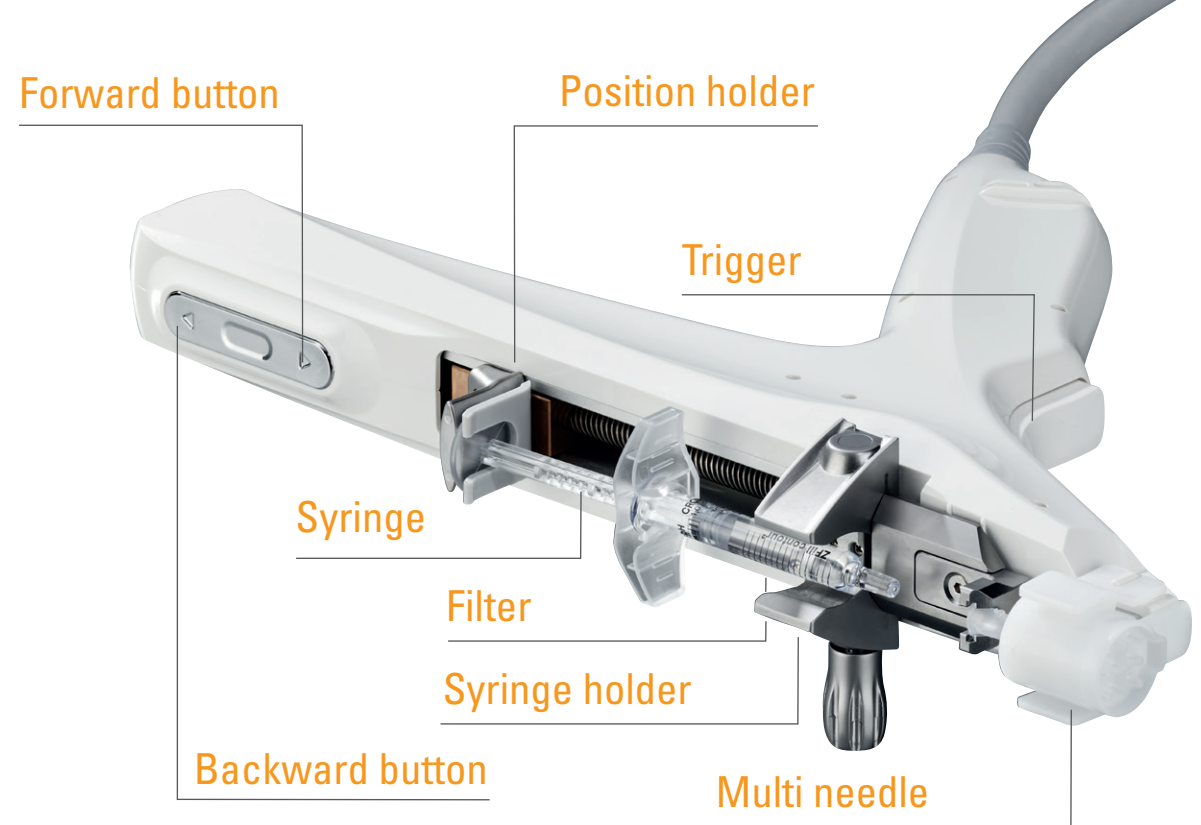

The concept

Our concept is combining a Filler treatment with our new reFill device that combines automatic injection and vacuum technology for precise injections and an overall uniform result.

Our reFill device can be adapted to different syringes and liquids.

FAST INJECTION

Automatic injection system

1

UNIFORM RESULT

Vacuum to ensure treatment depth

2

PRECISE TREATMENT

9-Pin needle

3

Combination

re|Fill

Better Compatibility, easier treatment

reFill has one handpiece and one multi needle adapter for a precise 9-pin injection.

All relevant parameters can be selected. Starting from the vacuum level, the syringe type and the specified dose to ensure a homogen result. The 9 pin needle adapter is single use to ensure the patient safety.





Homogenous result

Combines vacuum and automatic injection technology



Removable handpiece

Single use needle pin



Clear design for intuitive usage

- Different syringe types
- Selection of amount of injection

Scalp

Whole Face

The backs of the hands

Neck

Décollete

Upper knees

Application regions

one treatment - multiple effects



Perfect in form and function


Excellent Compatibility, easier treatment

7 inch LCD monitor with easy operation

Depth	Dose	Viscosity	Vacuum	HR mode
<div>DEPTH(mm) 2.0</div> <div>Detailed depth adjustment up to 5.0 mm by 0.2 mm interval</div>	<div>DOSE 10</div> <div>Detailed dose adjustment per one shot</div>	<div>VISCOSITY 3</div> <div>1~5 levels of speed control allow various solution treatments with different formulations</div>	<div>VACUUM 4</div> <div>1~5 levels of vacuum function apply to various body parts and customized treatment according to skin elasticity is available</div>	<div>MODE HR</div> <div>Vacuum function should be disabled for hair treatment.</div>



Operation principle of reFill

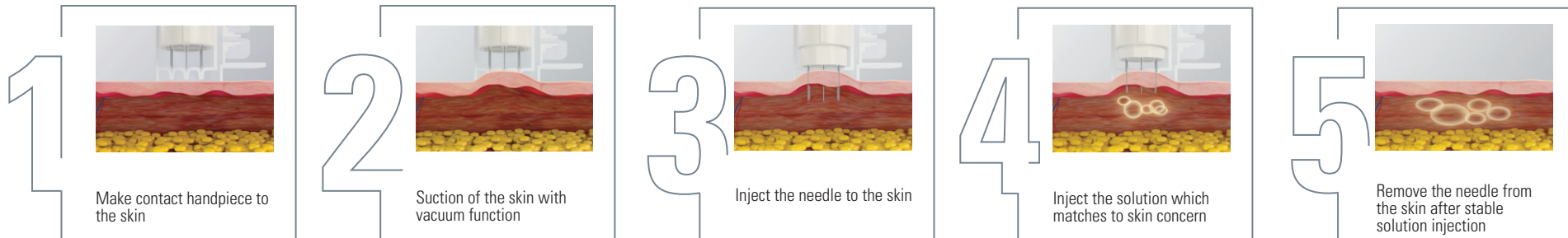

Technical data



**1, 2,
3, 5
CC**

Syringe Type

9pin needle

Needle Type

0.2 - 5.0 mm

Needle Depth

31G needle

Needle Gauge

AC 110-240 V, 50/60 Hz

Input Power

29 | 39 W (cm)

4 | 23 D (cm)

18 | 13 H (cm)

Dimensions Injector | Body

MD

0.9 | 4 kg

Weight Injector | Body



Follow us

Visit our homepage www.zimmer-aesthetics.com

